


SYSTEM CHART



*Add A or T to the article number for the desired obturator:
A220xxA standard obturator
A220xxT obturator with deflecting tip

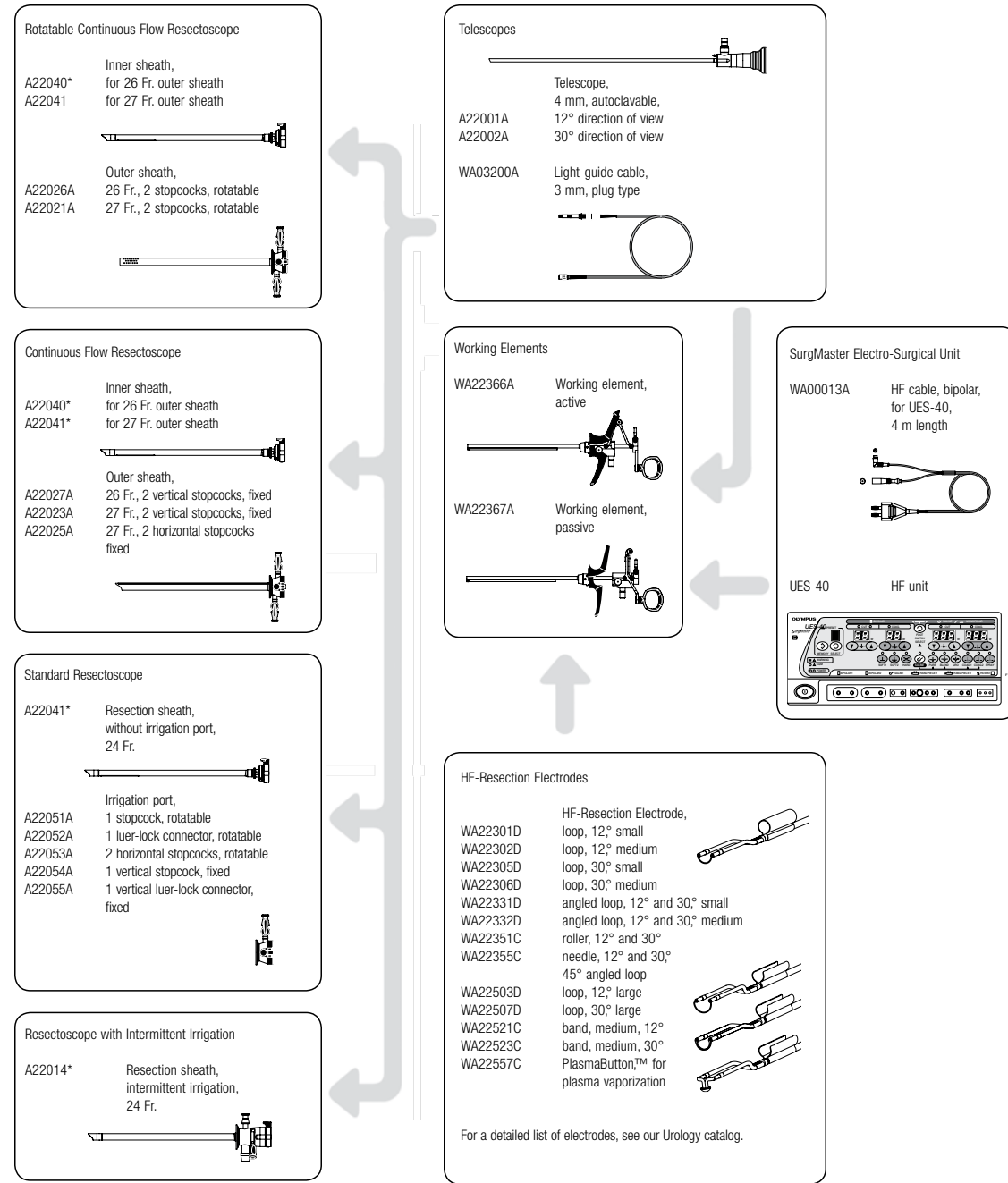

© 2010 Gyrus ACMI. All rights reserved. Printed in USA. TM and ® Trademarks and Registered Trademark, respectively, of Olympus and Gyrus ACMI, Inc. or their affiliated entities in the U.S. and/or other countries of the world. BR145 0910 Subject to change without notice.

For more information or order inquiries, please call 888.524.7266 visit us at www.gyrusacmi.com

Gyrus ACMI

136 Turnpike Road
Southborough, MA 01772
TEL 1-508-804-2600

Sales Locations:

USA

Gyrus ACMI

6655 Wedgwood Rd, Suite 160
Maple Grove, MN 55311
TEL 1-763-416-3000

Canada

Olympus Canada, Inc.
151 Telson Road
Markham, Ontario L3R 1E7
Canada
TEL 1-800-387-0437

Latin America

Olympus Latin America
5301 Blue Lagoon Drive, Suite 290
Miami, FL 33126, United States of America
TEL 1-305-266-2332

UK

KeyMed Ltd
KeyMed House
Stock Road
Southend-on-Sea, Essex SS2 5QH
United Kingdom
TEL 44-0-1702-616333

Europe

Olympus Europa Holding GmbH
Wendenstraße 14-18
20097 Hamburg
Germany
TEL 49-40-23773-0

Australia

Olympus Australia Pty Ltd.
31 Gilby Road
Mount Waverley
Victoria 3149, Australia
TEL 61-1300-132-992

Hong Kong

Olympus Hong Kong and China Ltd.
L43, Office Tower, Langham Place
8 Argyle Street, Mongkok
Kowloon, Hong Kong
TEL 85-2-2481-7812

Korea

Olympus Korea Co., Ltd.
4F, Gyeongang Bldg.
Kangnam-Gu, Seoul, 135-090
Korea
TEL 82-2-6207-3200

Singapore

Olympus Singapore PTE Ltd.
491B River Valley Road #12-01/04
Valley Point Office Tower
Singapore 248373
TEL 65-6834-0010

China


Olympus Sales & Service Co., Ltd.
R1202, NCI Tower, A12 Jianguomenwai Ave.
Chaoyang District
Beijing, 100022, China
TEL 010-65693535



Your Vision, Our Future

Plasma Vaporization Therapy

Resection in Saline with the PlasmaButton™ Electrode

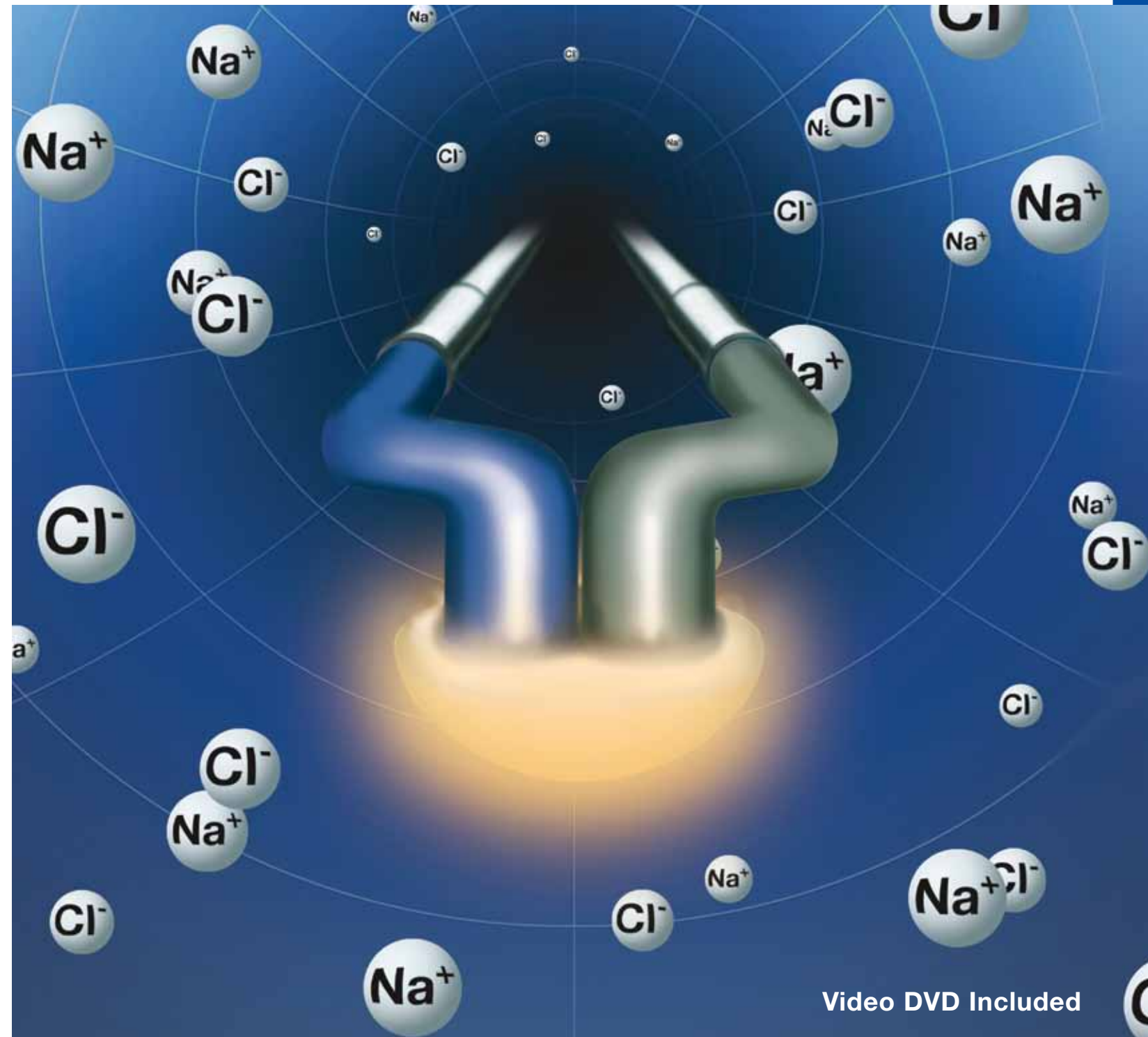

Video DVD Included

PLASMABUTTON™ ELECTRODE

The revolutionary, optimized spherical shape of the PlasmaButton™ Vaporization Electrode combined with the easy-to-learn vaporization technique results in well-coagulated, smooth tissue. Plasma Vaporization Therapy is safely performed in saline while the plasma corona reduces the penetration depth of the energy, affecting only localized tissue. The PlasmaButton™ electrode shape combined with resection-in-saline technology provides a safe, cost-efficient and easy-to-use solution to tissue management challenges.

- Virtually bloodless vaporization of the prostate
- Shorter procedure time
- No risk of hyponatremia
- Less bleeding (has even been demonstrated in patents on anti-coagulants).
- Short catheterization time
- Reduced risk of obturator nerve stimulation
- Minimal charring / thermal spread
- High ablation rate

WHAT IS PLASMA ENERGY?

Plasma is an electrically conductive gas cloud that, in this case, is created by the combination of radio frequency, energy and saline. It is made of vapor and charged particles called ions (atoms that have an electric charge). Due to its conductivity, the plasma allows the energy to cross at lower energy levels. This effect leads to lower operating temperatures and therefore less thermal spread.


Button electrode for plasma vaporization WA22557C



THE BIPOLAR OES PRO RESECTOSCOPE FOR RESECTION IN SALINE

The OES Pro resectoscope is designed with convenience and comfort in mind. The placement of key connections makes it extremely well balanced, and the finger spread on the working element is optimized to reduce fatigue. The high precision working element reduces friction force to a minimum and allows for a very smooth cutting experience. The logical locking system assures quick and secure assembly with one touch buttons for quick release.

These devices are not intended to be used in treating cancer of the prostate.

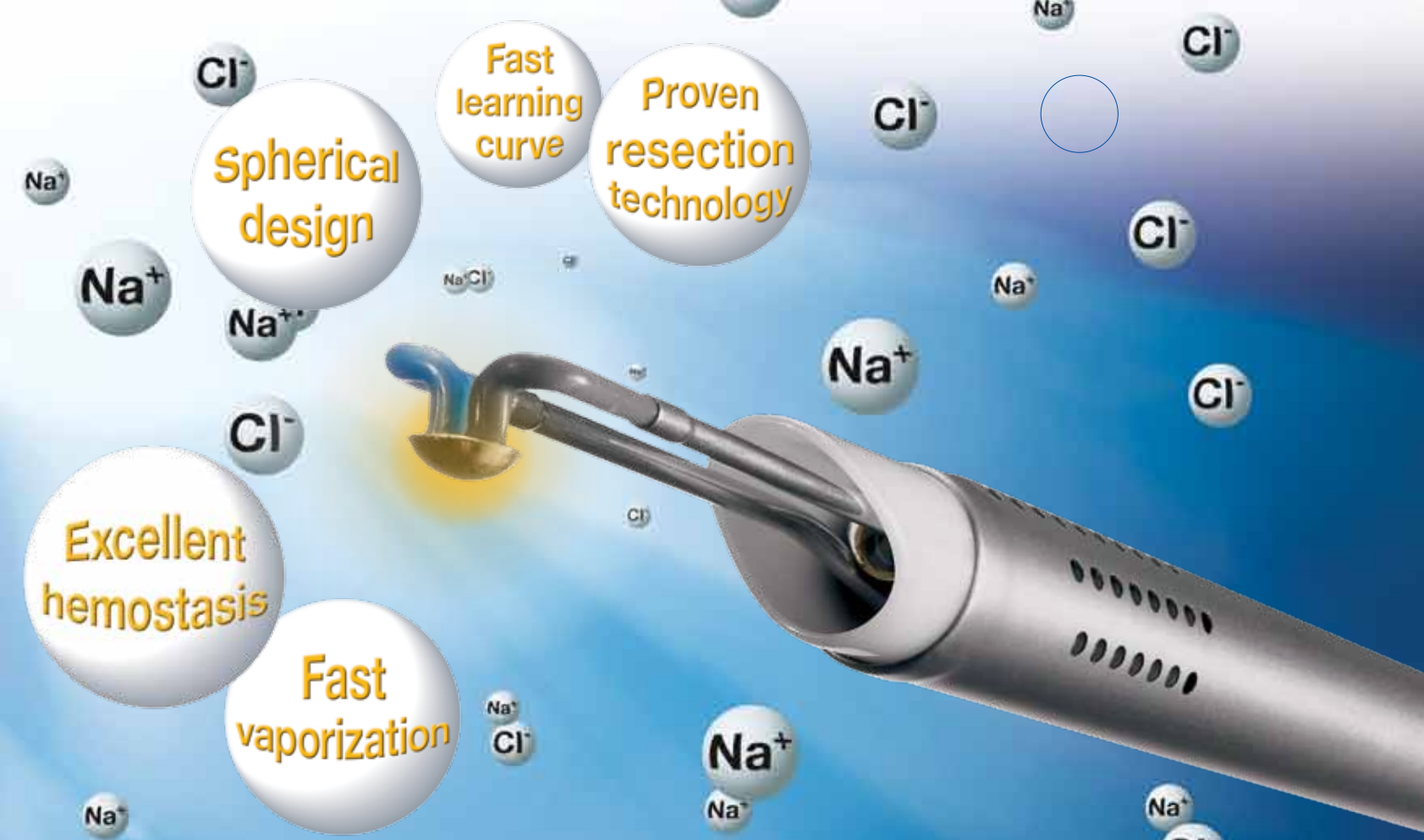

UES-40 ELECTROSURGICAL GENERATOR

The Olympus UES-40 SurgMaster — one generator for virtually any electrosurgical need. The Olympus SurgMaster covers almost every application where electrosurgery is performed. One unit enables the surgeon to perform resection in saline, conventional monopolar resection, and allows the connection of bipolar and monopolar hand instruments as well. SurgMaster — powered by Olympus.



- Resection in Saline Cable Outlet
- Pure Cut Mode — Minimum of 280 W
- Coag 1 Mode — Minimum of 140 W

UES-40 Generator

Software Version 2.0 or higher is required

